

BTS FREEEMG 300

Portable Surface EMG System
using Wireless Probes

BTS FREEEMG

High Performances Wireless 50 m - 9gr



BTS Bioengineering

Portable Surface EMG System using Wireless Probes



The most advanced Electromyography Unit in the world.

BTS FREEEMG 300 represents a generation leap in diagnostic device technology for surface dynamic electromyography analysis. Based entirely on wireless technologies, BTS FREEEMG 300 uses 16 miniaturized probes with active electrodes weighing less than 9 grams for signal acquisition and transmission.

BTS FREEEMG 300 stands for an improvement in the performance of the FREEEMG systems already appreciated by researchers and clinicians throughout the world. These practitioners need to be able to rely on a daily basis on a highly efficient system easy to configure and use, for differentiated clinical and researcher applications.

The probes amplify the EMG signals, digitize them and communicate with the compact and light receiving unit, which can be either way worn by the patient, left on the table or held by the doctor.

The complete absence of wiring not only minimizes the patient distress during the preparation but also grant him full range of motion during the task without restrictions.

The probes variable geometry and the substantially reduced size and weight allow them to be used on any body segment and during all types of movement (walking, running, jumping, etc) on a variety of subjects without affecting in any way the motor pattern.

Thanks to the use of environmentally friendly rechargeable batteries several hours of data can be recorded.

BTS EMG-Analyzer

BTS EMG-Analyzer is the most complete software solution for EMG signal analysis. Includes predefined templates for evaluations in clinic, sport and research field: jump, pliometry, gait, fatigue analysis, isokinetics, etc... It integrates also an editor to design elaboration protocols: the user can develop quickly and effectively customized elaboration protocols. Thanks to an innovative visual object interface, the mathematical analysis language is translated into graphical form.

Applications

BTS FREEEMG 300 can be used in research, sports, occupational medicine, gnathology, neurology and orthopaedics.

With the available software suites, BTS FREEEMG 300 becomes the most advanced diagnostic tool available today to evaluate:

- neurological and orthopaedic pathologies,
- pharmacological therapies,
- motor deficit progression,
- use of orthosis,
- rehabilitative follow-ups,
- sport training optimization

**Weigh less than 9 grams,
including battery**



US Partner:

zFlo, Inc
Boston, MA
508.269.5090



Real Dimension

Record of miniaturization

BTS FREEEMG 300 uses the best technology available today. Totally wireless, it integrates probes with active and variable geometry electrodes, weighing less than 9 grams.

Quicker Analysis and Increased Accuracy

The total absence of cables makes for a quick patient setup. The lightweight probes are attached directly to the pre-gelled electrodes and do not require any additional fixing such as adhesive tape. The patient can move around freely.

Powerful and Comprehensive

16 electromyography channels, 2 footswitch channels covering up to 8 contact areas, and up to 8 electrogoniometers or accelerometers probes.

16bit Accuracy

The sampling rate of up to 4 KHz and the 16 bit resolution give to the acquired signal the highest quality ever, providing "low noise" and the absence of movement artifacts.

On-board memory

The probes are equipped with a solid-state buffer memory to guarantee data safety in case of WiFi signal loss or when the probes are simply out of range during the acquisition.

More than 5h of non-stop acquisition

The large autonomy and energy-saver battery insure a whole day recording. Easy and quick recharging by clipping onto the proprietary battery charger. With "local" mode it can be used such as a holster, for continuously data monitoring and recording. For outdoor and out of the lab data collection applications.

Touch-screen interface

The system can be managed either from the workstation or from the mobile receiving unit, which is equipped with a high-resolution touch-screen display and an intuitive graphical interface.

Signal Range

Up to 50 metres for the transmission between the probes and the mobile receiving unit; up to 350 metres for the data transfer between the mobile unit and the workstation.


Interoperability

Designed to interoperate with movement analysis, posturography and stabilometry systems, thanks to the software integration through SDK or by means of the analogical output receiver. It is intended also to interoperate with isokinetic machines (BIODEX, CYBEX and CON-TREX), rehabilitative robots (Reo-go, BTS ANYMOV) and virtual reality therapeutic systems (BTS NIRVANA) for the evaluation of real muscle activity during rehabilitation or sport training and planned testing activities.

Quick setup



Signal range up to 350m

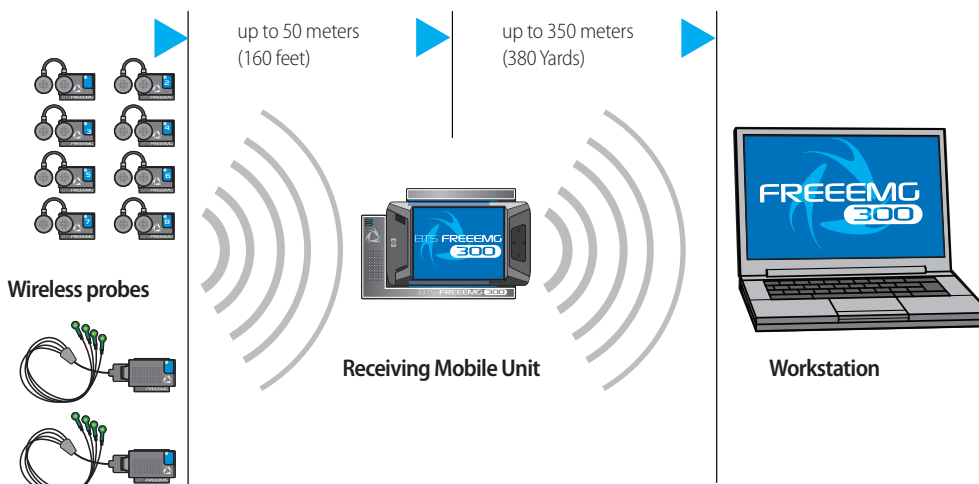

BTS FREEMG 300

Portable Surface EMG System using Wireless Probes



Technical features*

Wireless Probes

Surface electrodes	Variable geometry electrodes with mounting clip 16-bit resolution - acquisition frequency up to 4KHz
Data transmission	Wireless IEEE802.15.4 (probes - receiving unit)
Battery	Rechargeable with proprietary charger (clip connector)
Autonomy	5h battery life with non-stop recording, 5 days with Stand-by mode
Range probes - receiving unit	Up to 50 meters (160 feet) in free space (it can decrease to 30m -100 feet- in the presence of obstacles)
Memory	On board solid-state buffer memory system
Status LED	Acquisition/stand-by mode and low battery
Weight	<9 grams including battery
Size	23.8 x 37 x 10mm main electrode - Ø 16.5 x 10mm satellite electrode
Identification Labels	Available in four colors

Mobile Receiving Unit

EMG channels	Up to 16 wireless probes on each receiving unit
Data transmission	Wireless WiFi standard 802.11b (receiving unit – workstation)
Display	4"VGA touch-screen
Recording duration	Up to 9 hours with a single battery
Range receiving unit - workstation	Up to 30 meters (100 feet) indoor – up to 350 meters (380 yards) outdoors
Weight and size	360 grams – 155 x 106 x 35mm

Caution: U.S. Federal law restricts this device to sale by or on the order of a licensed healthcare practitioner

* Technical specifications are subject to change without prior notice.



Options

Software BTS EMG-Analyzer

Database for data storage. Graphic interface to build analysis template: operators for filtering, spectrum computation, contact event measure, latencies, thresholds, fix and mobile windows integration, RMS, interpolation, fatigue analysis, etc. Tool for the creation of customized report of the exams in HTML format completely configurable by the user.

BTS Docking Station

Multifunctional charger equipped with internal batteries to recharge the EMG probes when the power connection is not available.

BTS Workstation

Preconfigured Desktop or Laptop PC dedicated to the electromyography signal processing and to the wireless communications with the receiving unit.

Fine-Wire Adaptors

The EMG probes can be used also to acquire the EMG signal by means of Fine-Wire electrodes. The connection is easily realized through the dedicated adaptors equipped with clip and Fine-Wire links.

Wireless Foot Switch Probes

Independent FSR sensors for the automatic gait phases identification. Up to 4 single on-off sensors for each Foot Switch probe (8 in total).

Wireless Electrogoniometer Probes

Strain gauge technology for the accurate measurement of the angles drawn by the joints in the different planes.

Wireless Accelerometer Probes

Accelerometers for the Measurement of the linear acceleration in 3 axes. Integrated with Gyroscopes, permit kinematic measurement of position, velocity and acceleration. Allows the automatic gait phases identification (heel strike, toe contact) or the measurement of the amplitude and the frequency of a Tremor.

Video Acquisition System

BTS VIXTA, system to video record using up to 4 TV cameras simultaneously, natively synchronized with EMG signals

Analog Output Receiver

This tool allows for a quick integration with other motion analysis systems. It consists of a wireless receiver with analog output feature: raw EMG data available in analog format, simultaneously with the digital wireless data transmission. It is equipped with a USB interface for plug and play.

BTS FREEMG & ISOKINETIC Kit

For integration with isokinetic machines: BIODEX, CYBEX, CON-TREX.



BTS Bioengineering

HEADQUARTERS

VIALE FORLANINI 40
20024 GARBAGNATE MILANESE MI ITALY
TEL. +39 02.366.490.00
FAX +39 02.366.490.24

R&D CENTER

VIA DELLA CROCE ROSSA 11
35129 PADOVA PD ITALY
TEL. +39 049 981 5500
FAX +39 049 792 9260

BTS COMMERCIAL PARTNER NETWORK



US Partner:
zFlo, Inc
Boston, MA
508.269.5090